


CALZURO[®]

COLLECTION



CALZURO[®] Classic

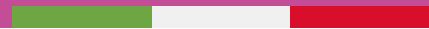

Professional clogs, made of a thermo-tecno polymer compound latex-free, completely antistatic, with side holes for foot perspiration, available with or without upper holes, autoclavable to 134°C/273°F, washable to 90°C/194°F. Slip-resistant sole CE certified. 100% Made in Italy.

Certified according to EN ISO 20347:2012, requirements:

- **A** (antistatic properties)
- **E** (energy absorption of seat region)
- **OB** (basic requirements)
- **SRA** (slip resistance)

MEDICAL DEVICE CLASS I NON STERILE: DM 1314685

SIZES

34/35 - 35/36 - 36/37 - 37/38 - 38/39 - 39/40 - 40/41 - 42/43 - 44/45 - 46/47

COLORS



WHITE



GREEN



LIGHTBLUE



ORANGE



BLACK



PINK ⁽²⁾



DARKBLUE



MET. BLUE



RED



PURPLE



MAROON



YELLOW



MET. GREY



TEAL ⁽¹⁾

(1) available only with upper holes

(2) available only up to 42/43

AVAILABLE ALSO WITH HEEL-STRAP



CECERTIFIED



ECO-FRIENDLY



100%RECYCLABLE



LATEX-FREE



ANTISLIP



AUTOCLAVABLE



ANTISTATIC



CRUELTY-FREE



MADEINITALY



WASHABLEIN90°
WASHINGMACHINE



UNISEX



IN ACCORDANCE
WITH REGULATION
(EU) 2016/425

ACCESSORIES Classic

INSOLE



Arch support insole for Calzuro, in special foam that allows feet to breathe, autoclavable to 121°C/250°F up to 50 times.

SIZES

34 - 35 - 36 - 37 - 38 - 39 - 40 - 42 - 44 - 46

CHARMS



CALZURELLI



SWAROVSKI



HEEL STRAP ON DEMAND



COLORS



WHITE



GREEN



LIGHTBLUE



ORANGE



BLACK



PINK



DARK BLUE



PURPLE



MAROON



YELLOW



MET. GREY



TEAL



RED



MET. BLUE

NB: FANCY COLORS SUPPLIED ONLY WITH WHITE HEEL STRAP

Autoclavable to 134°C/273°F

CALZURO[®] Fancy



Professional clogs, made of a thermo-tecno polymer compound latex-free, completely antistatic, with side holes for foot perspiration, available only with upper holes, autoclavable to 134°C/273°F, washable to 90°C/194°F.
Slip-resistant sole CE certified. 100% Made in Italy.

Certified according to EN ISO 20347:2012, requirements:

- **A** (antistatic properties)
- **E** (energy absorption of seat region)
- **OB** (basic requirements)
- **SRA** (slip resistance)

MEDICAL DEVICE CLASS I NON STERILE: DM 1314685

SIZES

34/35 - 35/36 - 36/37 - 37/38 - 38/39 - 39/40 - 40/41 - 42/43 - 44/45 - 46/47

COLORS



(1) available only with upper holes

(2) available only up to 42/43

SUPPLIED ONLY WITH WHITE HEEL-STRAP



CERTIFIED



ECO-FRIENDLY



100% RECYCLABLE



LATEX-FREE



ANTISLIP



AUTOCLAVABLE



ANTISTATIC



CRUELTY-FREE



MADE IN ITALY



WASHABLE IN 90°
WASHING MACHINE



UNISEX



IN ACCORDANCE
WITH REGULATION
(EU) 2016/425

CALZURO[®] Light



Professional clogs, made of ultra-light soft plastic latex-free material, casein-resistant, with antistatic plug and side holes for foot perspiration. Available only without upper holes and with back heel-strap. Machine washable up to 50°C/122°F. Slip-resistant sole CE certified. 100% Made in Italy.

Certified according to EN ISO 20347:2012, requirements:

- **A** (antistatic properties)
- **E** (energy absorption of seat region)
- **OB** (basic requirements)
- **SRC** (slip resistance)

MEDICAL DEVICE CLASS I NON STERILE: DM 1314687

SIZES

35 - 36 - 37 - 38 - 39 - 40 - 41 - 42 - 43 - 44 - 45 - 46

COLORS



WHITE



GREEN



NAVY BLUE



BLACK



PINK



MAUVE

NEW



LIGHTBLUE

NEW



OLIVE GREEN

(only up to size 41) (only up to size 44)



CERTIFICATICE



ECO-FRIENDLY



LATEX-FREE



ANTISCIVOLO



RESISTENTI
ALLACASEINA



ANTISTATICO



NONTESTATO
SUGLIANIMALI



MADENITALY



LA VABILIN
LAVATRICEA50°



UNISEX



IN ACCORDANCE
WITH REGULATION
(EU) 2016/425

ACCESSORIES **Light**



ANTISTATIC INSOLE

Comfortable ANTISTATIC insole, Drytech® fabric keeps the foot dry and latex foam bottom provides high level of comfort. Antibacterial treatment can prevent bacteria causing unpleasant odors.

SIZES

35 - 36 - 37 - 38 - 39 - 40
41 - 42 - 43 - 44 - 45 - 46



CUSTOM FIT INSOLE

Comfortable anti-perspiring insole, which **can be trimmed for custom fit**.

Drytech® fabric keeps the foot dry and latex foam bottom provides high level of comfort.

Antibacterial treatment can prevent bacteria causing unpleasant odors.

Machine washable up to 30°C/86°F



PILE INSOLE

Comfortable and warm PILE insole, which **can be trimmed for custom fit**.

Drytech® fabric keeps the foot dry and latex foam bottom provides high level of comfort. Antibacterial treatment can prevent bacteria causing unpleasant odors.

Machine washable up to 30°C/86°F

CALZURO[®] Buddy

Latex-free, antistatic slippers, washable up to 90°C/194°F, autoclavable to 121°C/250°F, with upper and side holes for foot perspiration, insulating sole from cold and hot, anatomical inner arch-support. Supplied with adjustable heel-strap. 100% Made in Italy.

SIZES

35 - 36 - 37 - 38 - 39 - 40/41 - 42/43 - 44/45

COLORS



WHITE +
METALLIC
PURPLE*

*only up to
size 42/43



WHITE +BLACK



WHITE +METALLIC
BLUE



100%RECYCLABLE



LATEX-FREE



AUTOCLAVABLE



ANTISTATIC



CRUELTY-FREE



MADEINITALY



WASHABLEIN90°
WASHINGMACHINE



UNISEX

CALZURO[®] ProSafety

(WITH SAFETY TOE-CAP)



Soft and practical professional footwear with closed back-heel, made of a brand-new special ultra-light copolymer, latex-free, casein-resistant.

Easy wear slip on and slip off and safe fit for maximum comfort and further your safety policy at work. Supplied without side vents, with antistatic plug and with synthetic polymer protection metal-free toe-cap (impact resistance of the toe-cap: 200 Joule on the toe).

The safety toe-cap provides adequate protection against mechanical risks, such as impact and crushing/breakthrough, but it shows greater lightness than the corresponding leather shoes with steel toe-cap.

Slip-resistant sole CE certified. 100% Made in Italy.

Certified according to EN ISO 20345:2011, requirements:

S4 (additional requirements):

- **Heel area closed**
- **A** (antistatic properties)
- **E** (energy absorption of seat region)
- **C leated outsole**

SRC (slip resistance)

MEDICAL DEVICE CLASS I NON STERILE: DM 1770379

SIZES

36 - 37 - 38 - 39 - 40 - 41 - 42 - 43 - 44 - 45 - 46

COLORS



WHITE



BLACK



CE CERTIFIED



ECO-FRIENDLY



ANTISLIP



LATEX-FREE



SAFETY TOE-CAP
METAL-FREE



CASEIN-RESISTANT



ANTISTATIC



CRUELTY-FREE



MADE IN ITALY



WASHABLE IN 50°
WASHING MACHINE



UNISEX



EASY-FIT



IN ACCORDANCE
WITH REGULATION
(EU) 2016/425

CALZURO[®] ProSP

(WITHOUT SAFETY TOE-CAP)



Soft and practical professional footwear with closed back-heel, made of a brand-new special ultra-light copolymer, latex-free, casein-resistant.

Easy wear slip on and slip off and safe fit for maximum comfort and further your safety policy at work.

Supplied with side vents for foot perspiration and with antistatic plug.

Slip-resistant sole CE certified. 100% Made in Italy.

Certified according to EN ISO 20347:2012, requirements:

- **A** (antistatic properties)
- **E** (energy absorption of seat region)
- **OB** (basic requirements)
- **SRC** (slip resistance)

MEDICAL DEVICE CLASS I NON STERILE: DM 1639834

SIZES

36 - 37 - 38 - 39 - 40 - 41 - 42 - 43 - 44 - 45 - 46

COLORS



WHITE



NAVY BLUE



RASPBERRY
(ONLY UP TO SIZE 42)



CE CERTIFIED



ECO-FRIENDLY



ANTISLIP



LATEX-FREE



CASEIN-RESISTANT



ANTISTATIC



CRUELTY-FREE



MADE IN ITALY



WASHABLE IN 50°
WASHING MACHINE



UNISEX



EASY-FIT



IN ACCORDANCE
WITH REGULATION
(EU) 2016/425

CALZURO[®] Bab

Antistatic white boots in special polymer compound (nitrile rubber), non toxic, resistant to hydrocarbons (fuel oil). A good washing resistance at medium-high temperatures 85°C/185°F helps maintain within optimal limits to the spread of bacteria. Resistant to animal fats, vegetable oils and casein, slip-resistant sole. Available: model at knee.

White polyester knitting inside.

Washable up to 85°C/185°F. 100% Made in Italy.

BAB (to the knee)

Certified according to EN ISO 20347:2012, requirements:

- O4 (closed seat region, antistatic properties, energy absorption of seat region)
- FO (bottom sole resistance to fuel oil)
- SRC (Slip resistance)

SIZES

36 - 37 - 38 - 39 - 40 - 41 - 42
43 - 44 - 45 - 46 - 47 - 48

COLORS



WHITE



CERTIFIED



ANTISLIP



FUELOIL-RESISTANT



CASEIN-RESISTANT



ANTISTATIC



CRUELTY-FREE



MADE IN ITALY



UNISEX



WASHABLE IN 85°
WASHING MACHINE



White boots in special polymer compound (nitrile rubber), non toxic, model at calf, white polyester kitting inside.
Washable up to 85°C/185°F. 100% Made in Italy.



BAB Half Wellington (to the calf)

Certified according to EN ISO 20347:2012, requirements:

- **OB** (basic requirements)
- **SRA** (slip resistance)

SIZES

36 - 37 - 38 - 39 - 40 - 41
42 - 43 - 44 - 45 - 46 - 47

COLORI



WHITE



CE CERTIFIED



MADE IN ITALY



UNISEX



ANTISLIP



FUEL OIL-RESISTANT



CRUELTY-FREE



WASHABLE IN 85°
WASHING MACHINE

CALZURO[®] Effetto

Ciabatte per il personale medico e paramedico, suola in PU (poliuretano) antistatico ed antiscivolo certificata, tomaia in pelle pigmentata e perforata per una traspirazione sana e naturale del piede, morbida imbottitura sul collo del piede. Possono essere fornite con cinturino posteriore.

Certified according to EN ISO 20347:2012, requirements:

- **A** (antistatic properties)
- **E** (energy absorption of seat region)
- **OB** (basic requirements)
- **SRC** (slip resistance)

MODELS



WITHOUT HEEL-STRAP



WITH HEEL-STRAP



SIZES

35 - 36 - 37 - 38 - 39 - 40 - 41
42 - 43 - 44 - 45 - 46

COLORS



WHITE
(white sole)



BLUE
(black sole)



ANTISTATIC



UNISEX



MADE IN ITALY



CE CERTIFIED



ANTISLIP

CALZURO[®] Giulia

Women's leather slippers with removable double-face insole* (massaging/leather) and upper heel-strap to adjust the instep, TR (thermoplastic material) bottom sole.

MODELS



WITH UPPER HOLES
AND WITHOUT HEEL-STRAP



WITH UPPER HOLES
AND WITHOUT HEEL-STRAP



WITHOUT UPPER HOLES
AND WITH HEEL-STRAP



*ALWAYS SUPPLIED WITH
FOOTWEAR ONLY

SIZES

35 - 36 - 37 - 38 - 39 - 40 - 41

COLORS



WHITE
(white sole)



BLUE
(Available only with upper holes
and without heel-strap)



CALZURO[®] Wood

Calzuro[®] Wood the “flexible”, with perforated upper in pure Italian leather and anti-bacterial fabric inside, flexible sole in beech wood and special rubber. Washable up to 40°C/104°F.

SIZES

35 - 36 - 37 - 38 - 39 - 40
41 - 42 - 43 - 44 - 45 - 46

COLORS



WHITE



NAVY BLUE



WASHABLEIN40°
WASHINGMACHINE



UNISEX

CALZURO[®] Comfort

Comfort Slippers entirely made of leather and sewed by hand, with back heel-strap, anatomical insole in natural perspiring leather, TR (thermoplastic material) bottom sole.

SIZES

35 - 36 - 37 - 38 - 39 - 40 - 41

COLORS



WHITE



BLUE



CALZURO[®] Aqua

Flip flops made of EVARITE[®], a new, soft, lightweight, latex-free plastic material which does not favor the proliferation of fungi and bacteria.

Anatomically shaped footbed, well-dressing anatomic shape, thong comfortable to wear, slightly raised heel and toe bed.

Machine washable up to 50°C/122°F. Slip-resistant sole CE certified. 100% Made in Italy.

Certified according to EN ISO 20347:2012, requirements:

- E (energy absorption of seat region)
- OB (basic requirements)
- SRC (slip resistance)

MEDICAL DEVICE CLASS I NON STERILE: DM 1587673

SIZES

35 - 36 - 37 - 38 - 39 - 40
41 - 42 - 43 - 44 - 45 - 46

COLORS



CE CERTIFIED



ECO-FRIENDLY



LATEX-FREE



ANTISLIP



CRUELTY-FREE



MADE IN ITALY



WASHABLE IN 50°
WASHING MACHINE



UNISEX

HOW TO SANITIZE CALZURO® CLOGS

“CLASSIC” TYPE



MODEL “S” (WITH OR WITHOUT UPPER HOLES)

- Sterilizable in autoclave up to 134°C / 237°F
- Washable in clogwasher up to 90°C / 194°F with specific disinfecting cleaner
- For immersion in 5% sodium hypochlorite aqueous solution for 15 minutes and subsequent rinsing
- By spraying of pure alcohol at 95% and leaving it for at least 15 minutes (optional subsequent rinsing)
- For Immersion with specific disinfectant (such as: Napisan®), thoroughly following the manufacturer’s instruction and timing
- Hand washing at 30°C minimum for at least 3 minutes with whichever soap and leaving it in for at least 15 minutes before rinsing



MOD.BASE (SENZA FORI SUPERIORI)

- sterilizzabile in autoclave fino a 121°C
- lavabile in lavazoccoli fino a 90°C con detergenti disinfettanti dedicati
- per immersione in soluzione acquosa al 5% di ipoclorito di sodio per 15 minuti e successivo risciacquo
- per spruzzamento di alcol puro a 95° lasciandolo agire per almeno 15 minuti (successivo risciacquo è opzionale)
- in immersione con disinfettanti dedicati (esempio: Napisan®) con istruzioni e tempi secondo dati fabbricante
- con lavaggio a mano ad almeno 30°C per almeno 3 minuti con qualsiasi sapone, lasciandolo agire per almeno 15 minuti prima del risciacquo

MODELS IN EXPANDED EVA



(CALURO LIGHT – CALZURO PRO – CALZURO PRO SAFETY – AQUA)

- Washable in clogwasher up to 50°C / 122°F with specific disinfecting cleaner
- For immersion in 5% sodium hypochlorite aqueous solution for 15 minutes and subsequent rinsing
- By spraying of pure alcohol at 95% and leaving it for at least 15 minutes (optional subsequent rinsing)
- For Immersion with specific disinfectant (such as: Napisan), thoroughly following the manufacturer’s instruction and timing
- Hand washing at 30°C minimum for at least 3 minutes with whichever soap and leaving it in for at least 15 minutes before rinsing





100% MADE IN ITALY

www.calzuro.it

BIHOS SAS DI ROBERTA BIDOIA
Via Dell'Artigianato, 18
35010 VIGONZA (PD)

C.F. e P.IVA 01865630287
Tel. 049-8930259

www.bihos.com
bihos@bihos.com
info@calzuro.it

Pat Coakley
Medical 

Pat Coakley Medical Ltd.

Units E3 E4 Network Enterprise Park
Kilcoole A63 H244 Co Wicklow
IRELAND

+353 12018070

sales@patcoakleymedical.ie