




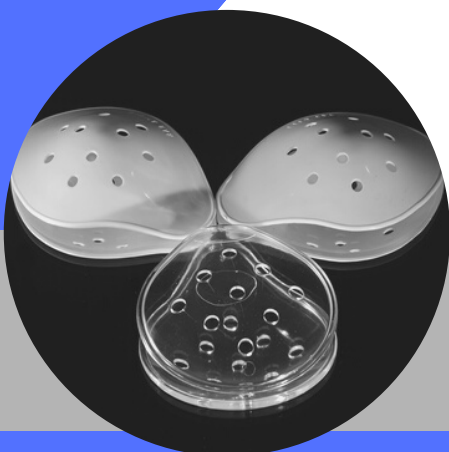


INSIGHT
Medical Products Ltd

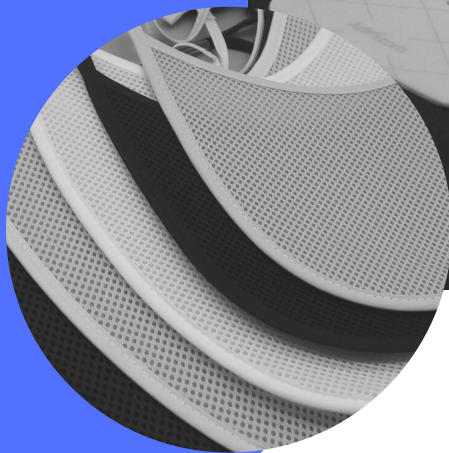

PRODUCT CATALOGUE

2024

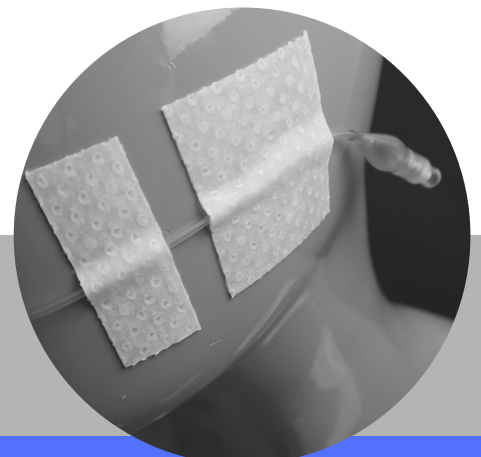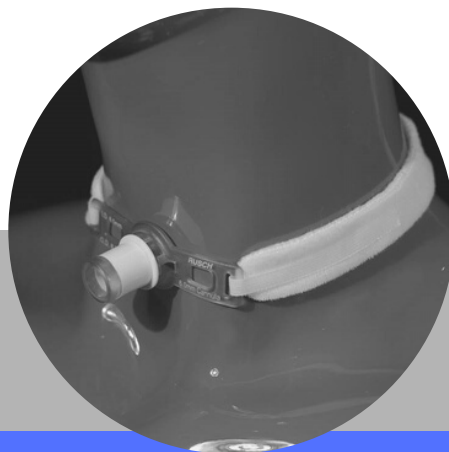

AN INSIGHT INTO INSIGHT

Insight Medical Products Ltd. was established by its current directors in May 2000. From our base in the Cotswolds, Insight Medical has become a highly respected supplier of good value, high quality innovative products to the healthcare market.

We are pleased to work with healthcare professionals to develop and manufacture innovative products that will help budgets and improve patients well-being.

At Insight Medical we have our own in-house manufacturing facility with a fully accredited cleanroom facility for manufacturing and packing.



Our non-sterile medical devices comply to the Medical Device Regulations (EU) 2017/745 and applicable harmonised standards. Our sterile devices meet the requirements to the Medical Device Directive 93/42/EEC, UK MDR 2002 and applicable harmonised standards under supervision of SGS. We are currently in transition for the sterile devices to be MDR Certified.

Our Quality Management System complies and meets the requirements of ISO 13485:2016 under supervision of SGS.

We are in progress of obtaining ISO 14001 -Environmental Management Systems certification under supervision of BAB.

We pride ourselves on our reputation for Customer Service and are committed to ensuring a first class service to all of our customers, in the United Kingdom and overseas.



TRACHEOSTOMY CARE

- Tracheostomy Tube Holder
- Tracheostomy Speaking Valve
- Tracheostomy Tube Cleaning Swab
- Advazorb Plus (T) Tracheostomy Dressing
- Advazorb Lite (T) Tracheostomy Dressing
- Advadraw (T) Laminated Tracheostomy Dressing
- Advazorb Fixation Tracheostomy Dressing

TUBE FIXATION

- Endotracheal Tube Holder
- Insil Silicone Tape
- Transducer Holder
- Tube Tidy Clips

LARYNGECTOMY PROTECTORS

- Sofshield Laryngectomy Protector
- Cascade Shower Protector

LARYNGECTOMY CARE

- Purifoam Laryngectomy Filter
- Voice Prosthesis Brush

THEATRE CONSUMABLES

- Eye Shields
- Bite Guard

MISCELLANEOUS

- Amniotic Perforator

VETERINARY

- Cat Bladder Irrigation Catheter
- Dog Bladder Irrigation Catheter

TRACHEOSTOMY CARE

TRACHEOSTOMY TUBE HOLDER



Our Tracheostomy Tube Holder uses high quality materials to provide a safe and secure fixation for the tracheostomy tube.

The neckband is manufactured from an ultra soft brushed cotton/foam material for prevention of soreness and discomfort to the patient.

A simple to use two-piece design allows quick fitting to the majority of neck sizes, while the integral 'cough strip' on the adult variant enables reflex movement to the patient if required.



- Adult and Paediatric sizes available
- 2-Piece design allows secure fitting for most devices
- Cough reflex section offers more patient comfort
- Available on FP10 Prescription
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

The Tracheostomy Tube Holder is used to secure the tracheostomy tube and to aid prevention of accidental decannulation for patients with a tracheostomy tube in place.

Safety Information

- The holders cannot be used on silicone tracheostomy tubes
- Change the holder daily or more frequently if required
- Discard the device should there be signs of deterioration
- The holder must be replaced each time the tracheostomy tube is cleaned
- The holder must be changed each time after patient bathing
- Single patient use only

Medical Symbols



LATEX FREE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION

NHS
Supply Chain

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
TH/100	Tracheostomy Tube Holder	Adult	289-4814	FTH000	10
TH/150	Tracheostomy Tube Holder	Paediatric	289-4806	FTH092	10

TRACHEOSTOMY CARE

TRACHEOSTOMY SPEAKING VALVE




The Trachie-Voice low profile Tracheostomy Speaking Valve fits securely to a 15mm connector tracheostomy tube.


The device has a one way mechanism that opens when the patient breathes in, allowing air into the airway via the tracheostomy.

The valve closes on expiration allowing air to follow the normal route of expiration through the mouth and nose, instead of the tracheostomy, permitting speech.

Being 'hands free' allows voicing to take place without the need for finger occlusion, so reducing possible infection and hygiene concerns.

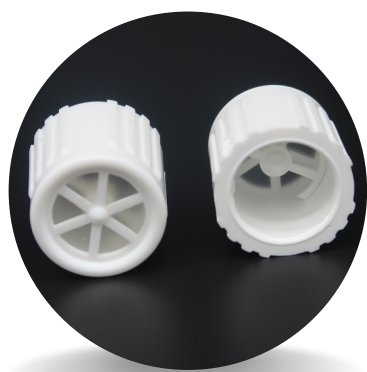

- Available in two colour variants
- For non-ventilated patients
- Economical/low cost device
- Lightweight to reduce pressure on the stoma site
- Reduces hygiene and infection concerns
- Ridged sides enabling quick application and removal
- Supplied complete with storage box
- Can be used for up to 30 days
- Available on FP10 Prescription
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

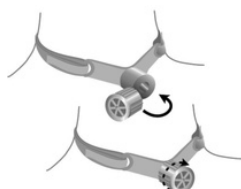

Intended Purpose

The Tracheostomy Speaking Valve is developed to help with the treatment of Tracheostomy patients. The valve provides patients the ability for voicing and to communicate more effectively. Speaking valves should never be used with the tracheostomy tube cuff inflated as this can lead to asphyxiation and respiratory arrest.

Safety Information

- The tracheostomy cuff must always be deflated when in use.
- If the flexible valve does not open during inspiration or close during expiration, remove the Speaking Valve immediately.
- The initial fitting of the Speaking Valve should be carried out by a trained clinician or caregiver and observed for a period to ensure adequate respiration.
- The device is not for use with patients who are being mechanically ventilated.
- Discard the device should there be signs of deterioration.
- Single patient use only

Medical Symbols



LATEX FREE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION

NHS
Supply Chain


IMP CODE	DESCRIPTION	COLOUR	DRUG TARIFF	NHS CODE	PACK QTY
TSV/100	Tracheostomy Speaking Valve	White	410-7546	FDG1636	12 (6x2)
TSV/150	Tracheostomy Speaking Valve	Blue	-	FDH721	12 (6x2)

TRACHEOSTOMY TUBE CLEANING SWAB

Our Tracheostomy Cleaning Swabs are intended for cleaning all types of PVC and silicone inner-tubes from tracheostomy or laryngectomy tubes.

Available in two head sizes, our Tube Cleaning Swabs are secure and suitable for popular sizes of tubes up to 6mm (small) and 9mm (large and extra long).

The foam heads of the swabs are heat welded for added security whilst cleaning. The soft foam heads help reduce damaging or scuffing the inside of the tracheostomy tube whilst cleaning. The white plastic handle is flexible for ease of use during cleaning.

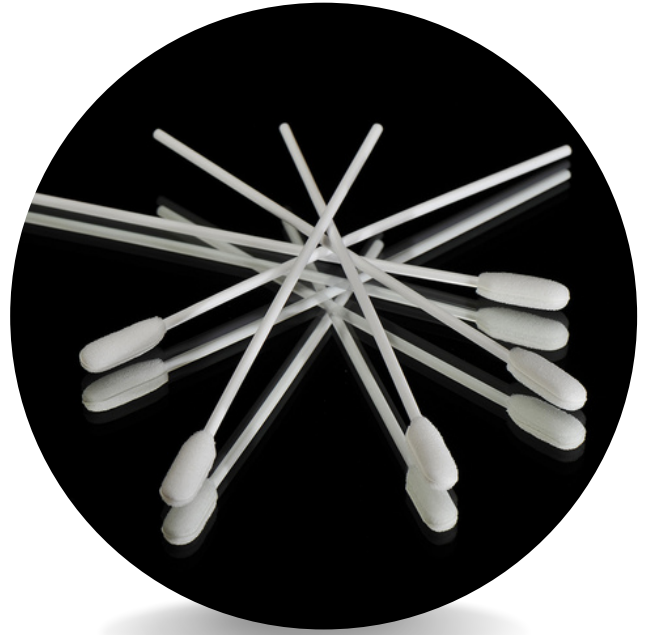

- Use with warm water to remove mucus residue in the tube
- Sizes available for popular sizes of tube
- Supplied in a resealable tray
- Available on FP10 Prescription
- Latex Free
- CE Marked

Intended Purpose

The Tracheostomy Tube Cleaning Swabs are used for cleaning all types of inner cannula tracheostomy tubes. The swabs are developed to prevent damage to the tubes which could lead to contamination.

Safety Information

- Inner cannula tracheostomy tubes must be cleaned whilst not in situ
- Refer to the tracheostomy tube manufacturer's instructions before cleaning
- Always ensure the inner cannula is clean before use, before inserting back into the patient's tracheostomy
- Do not excessively bend the handle of the swab
- Only for the intended use on inner cannula tracheostomy tubes and not for use on other devices
- Single patient use only

Medical Symbols



LATEX FREE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION

NHS
Supply Chain

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
TCS/900	Tracheostomy Tube Cleaning Swab	Large	324-8689	FDH420	100 (10x10)
TCS/950	Tracheostomy Tube Cleaning Swab	Small	324-8671	FDH421	100 (10x10)
TCS/1000	Tracheostomy Tube Cleaning Swab	Extra Long	-	FDH611	100 (10x10)

IMP CODE	HEAD WIDTH	HEAD LENGTH	HANDLE LENGTH	MAXIMUM TUBE SIZE
TCS/900	13mm	26mm	155mm	9mm
TCS/950	8mm	22mm	155mm	6mm
TCS/1000	9mm	25mm	245mm	9mm

TRACHEOSTOMY CARE

ADVAZORB PLUS (T) TRACHEOSTOMY DRESSING





Advazorb Plus (T) offers a good standard in hydrophilic foam tracheostomy dressings with high absorbency ability.

Manufactured in lightweight soft, and conforming low-adherent foam with high fluid absorbance and retention characteristics.

Manufactured with a film backing to ensure improved strength and absorbance for medium/high exuding wounds.

Supplied sterile and individually wrapped, ready for use. Available in one convenient size to suit most intubation procedures, and to offer comfort to the patient.

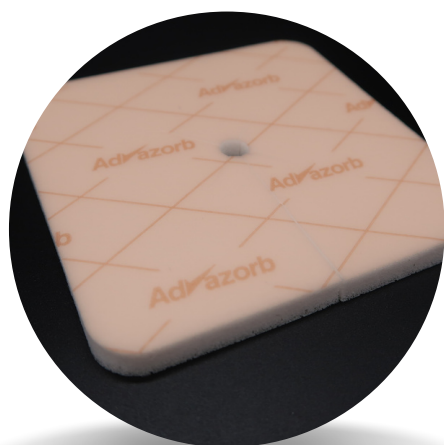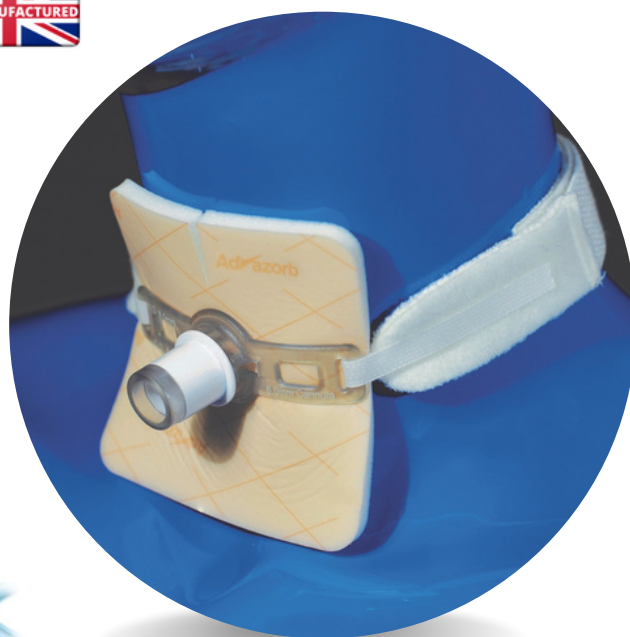

- Cushioning effect on the wound
- Low adherence -Reduced trauma on removal
- Highly absorbent 5mm foam
- Excellent patient comfort
- Rapid uptake of exudates
- Available on FP10 Prescription
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

Advazorb Plus Tracheostomy Dressing is developed to help with the treatment for Tracheostomy patients.

Safety Information

- Do not use with oxidising agents such as hypochlorite solutions or hydrogen peroxide as these can break down the absorbent component of the dressing
- Not suitable for arterial bleeds or heavy bleeding wounds
- Do not use if packaging is damaged as sterility may be reduced
- Single patient use only

Medical Symbols



STERILE

LATEX FREE



Distributed by Insight Medical Products Ltd.

ORDERING INFORMATION

NHS
Supply Chain

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/3900	Advazorb Plus (T) Dressing	80x80mm	395-1126	FDH419	20

TRACHEOSTOMY CARE

ADVAZORB LITE (T) TRACHEOSTOMY DRESSING





Advazorb Lite (T) is the definitive lightweight hydrophilic slim foam tracheostomy dressing.

The tracheostomy dressing offers remarkable standards of conformability, softness and lightness compared to gauze dressings.

Intended for the treatment of light/medium exuding wounds as moisture is absorbed and retained in the hydrophilic cellular structure.

Supplied sterile and individually wrapped, ready for use. Available in one convenient size to suit most intubation procedures, and to offer comfort to the patient.

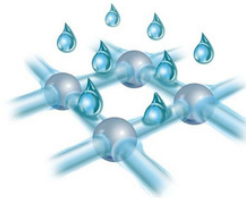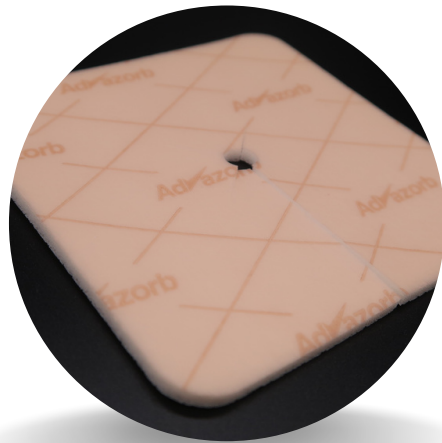

- Cushioning effect on the wound
- Low adherence -Reduced trauma on removal
- Highly absorbent 2.5mm foam
- Excellent patient comfort
- Rapid uptake of exudates
- Available on FP10 Prescription
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

Advazorb Lite Tracheostomy Dressing is developed to help with the treatment for Tracheostomy patients.

Safety Information

- Do not use with oxidising agents such as hypochlorite solutions or hydrogen peroxide as these can break down the absorbent component of the dressing
- Not suitable for arterial bleeds or heavy bleeding wounds
- Do not use if packaging is damaged as sterility may be reduced
- Single patient use only

Medical Symbols



STERILE

LATEX FREE



Distributed by Insight Medical Products Ltd.

ORDERING INFORMATION

NHS
Supply Chain

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/4212	Advazorb Lite (T) Dressing	80x80mm	374-0230	ELY392	20

TRACHEOSTOMY CARE

ADVADRAW (T) LAMINATED TRACHEOSTOMY DRESSING



Advadraw (T) is a tri-laminate tracheostomy dressing which has a low adherent perforated wound film contact layer, an absorbent centre action and a strike-through resistant layer to seal the dressing.

The tracheostomy dressing is comprised of an enhanced lighter fibrous material, with tests showing improved absorbency rates when compared to other manufacturers.

The slim contours of Advadraw (T) avoids elevating the tracheostomy tube away from the neck, while the curved edges ensure maximum comfort to the patient.

Supplied sterile and individually wrapped, ready for use. Available in one convenient size to suit most intubation procedures, and to offer comfort to the patient.



- Cushioning effect on the wound
- Low adherence -Reduced trauma on removal
- Conformable patient protection
- Suitable for low/medium levels of exudates
- Rapid uptake of exudates
- Available on FP10 Prescription
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

Advadraw Tracheostomy Dressing is developed to help with the treatment for Tracheostomy patients.

Safety Information

- Advadraw should be used with caution in the treatment of a highly exuding stoma or leaking drainage site
- Not suitable for arterial bleeds or heavy bleeding wounds
- Do not use if packaging is damaged as sterility may be reduced
- Single patient use only

Medical Symbols



STERILE

LATEX FREE



Distributed by Insight Medical Products Ltd.

ORDERING INFORMATION

NHS
Supply Chain

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/4416	Advadraw (T) Dressing	80x80mm	409-4970	FDH615	20

TRACHEOSTOMY CARE

ADVAZORB FIXATION TRACHEOSTOMY DRESSING





Advazorb Fixation is a post tracheostomy decannulation dressing fitted with a central location button that prevents unwanted escape of air through the wound site with gentle pressure applied during occasions such as talking, coughing and sneezing.

The device comprises of a hydrophilic foam dressing with a unique perforated silicone adhesive wound contact layer. It is an atraumatic foam dressing which does not adhere to moist wounds.

The silicone helps maintain a moist wound healing environment and makes it ideal as an adherent wound contact dressing. Soft silicone dressings do not damage newly formed granulating tissue and do not lose their adhesion after initial application.

The silicone layer prevents the need for additional dressing or fixation. Silicone adhesive combines security with pain-free dressing changes. The waterproof backing enables the patient to shower without the need to replace or cover the dressing.

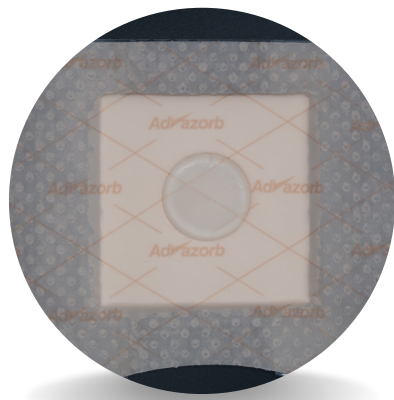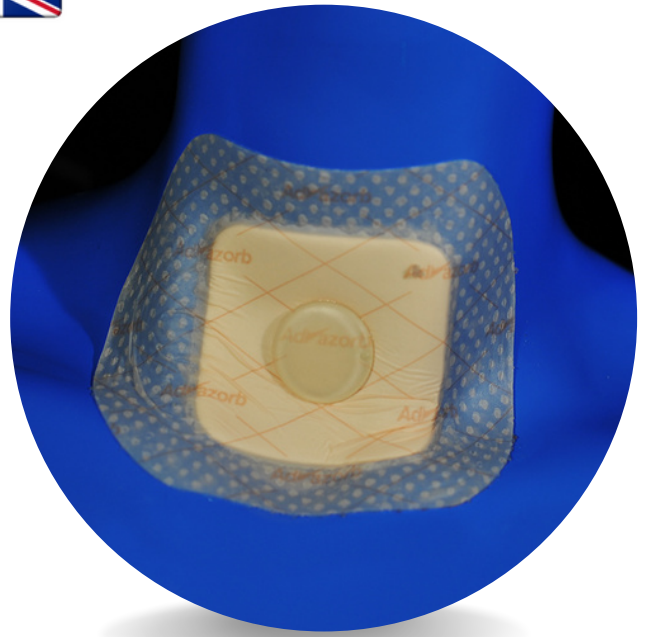

- Supplied sterile and individually wrapped, ready for use
- Location button enables accurate identification of the tracheostomy site
- Maintains a moist environment and manages exudates
- Allows atraumatic dressing change
- Provides a barrier for bacteria
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

Advazorb Fixation Tracheostomy Dressing is developed to help with the treatment for Tracheostomy patients.

Safety Information

- Do not use if allergic to silicone
- Cannot be used for arterial bleeds or heavy bleeding wounds
- Do not use if packaging is damaged as sterility may be reduced
- Single patient use only

Medical Symbols



STERILE

LATEX FREE



Distributed by Insight Medical Products Ltd.

ORDERING INFORMATION

NHS
Supply Chain

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/4341	Advazorb Fixation Dressing	12x10cm	-	FDH605	10

TUBE FIXATION

ENDOTRACHEAL TUBE HOLDER



The Endotracheal Tube Holder has been commended by professionals on the security and patient comfort of the device.

Manufactured with high quality soft materials, the product has proven to greatly reduce incidents of sores to the patients mouth, while maintaining a secure fitting to the ultra soft neckpiece.

The neckpiece is fitted with an adjustable laminated hook/loop tab which can be repositioned to size.

Each set is supplied with 3 securing straps to enable the user to replace soiled items, therefore improving value over alternative methods.



- Fully adjustable to fit most neck sizes
- Two sizes available
- Spare white neckbands available to order
- Laminated hook/loop tabs offer security of fitting
- Ultra soft materials reduce abrasions to the patient
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

The Endotracheal Tube Holder has been developed to secure the endotracheal tube when in use on patients during endotracheal intubation treatment, whilst also preventing case of patient mouth sores.

Safety Information

- Before use, clean tube with an alcohol wipe. Not providone-iodine
- Replace the blue securing strap when too wet or soiled with a spare strap provided
- Please refer to the training poster for guidance on using the device.
- Discard the device should there be signs of deterioration
- Single patient use only

DETAILED TRAINING POSTER AND PRODUCT TRAINING AVAILABLE ON REQUEST

Medical Symbols



LATEX FREE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION

NHS
Supply Chain

IMP CODE	DESCRIPTION	SIZE	NHS CODE	PACK QTY
TFP/200	Endotracheal Tube Holder Standard	20" Neckband	FTH005	10
TFP/201	Endotracheal Tube Holder Extra Long	24" Neckband	FTH025	10
TFP/222	Endotracheal Tube Holder (Neckband Only)	20" Neckband	FTH130	10
TFP/225	Endotracheal Tube Holder Extra Long (Neckband Only)	24" Neckband	FTH131	10

TUBE FIXATION

INSIL SILICONE TAPE



Insil is breathable tape which is backed with a silicone adhesive. When the tape is removed from application, it is extremely gentle on the skin. The tape can be repositioned if required.

The silicone tape can be used to fix most medical devices and for tube fixation and theatres for securing closure of eyelids during operations, and allowing painless removal.

Insil is available in 2cm and 4cm widths with fine perforations every 5cm to allow separation without the need of scissors. Supplied in pocket-sized dispenser boxes for ease of use



- Gentle adhesive qualities upon removal
- Available on FP10 Prescription
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

Insil Silicone Tape is primarily used as fixation tape for most medical devices or tubes. The device has a breathable soft side which is backed with silicone adhesive that enables gentle adhesive qualities upon removal of the tape. The upper side of the tape is soft for patients, providing comfort.

Safety Information

- Do not use if allergic to silicone
- The device is only intend for fixation of most lightweight medical devices and tubes. Insil Silicone Tape is not intended for wound care.
- Single patient use only

Medical Symbols



LATEX FREE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	PACK QTY
INSIL/200	Insil Silicone Tape	20mm x 3m	342-4132	12 (12x1)
INSIL/400	Insil Silicone Tape	40mm x 1.5m	342-4140	12 (12x1)

TUBE FIXATION



TRANSDUCER HOLDER



A simple and effective method for ensuring that the transducer remains in the correct position either within the department or while the patient is being moved and replaces the need for adhesive tape application.

The holder ensures stable positioning of the transducer for extended periods of time whilst providing patient comfort.

Made from soft material the transducer is attached securely to the holder by using an additional hook/loop strip through the eyelets of the flange of the transducer. The holder reduces the difficulties in taping the transducer to the patient or repositioning the transfusion stand due to patient movement, also saving nursing time.

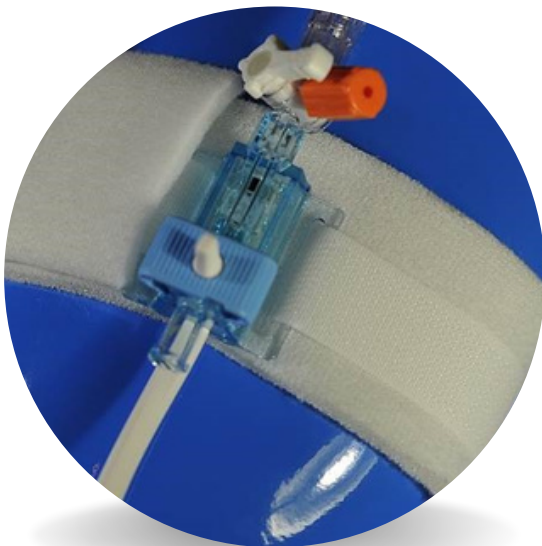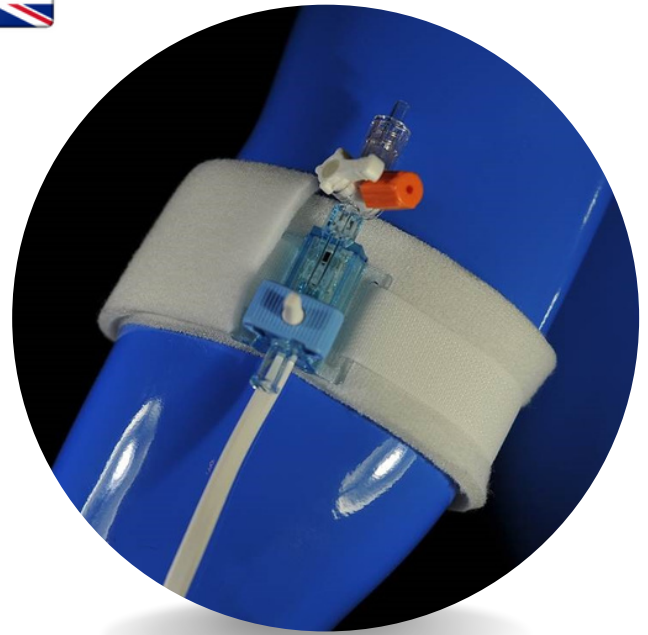

- For use in the department or during patient transfer
- Eliminates the use of adhesive tape on the patient
- Quality soft material for patient comfort
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

The Transducer Holder is used for securing Transducers in place on patients undergoing hemodynamic monitoring. The Transducer Holder when placed correctly by medical professionals, ensures the accurate level of the transducer at the correct phlebostatic axis to aid prevention for inaccurate blood pressure readings from the Transducer.

Safety Information

- Check the Transducer fits correctly before use
- Discard the device should there be signs of deterioration
- Single patient use only

Medical Symbols



LATEX FREE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	PACK QTY
TFP/400	Transducer Holder	Standard	-	10

TUBE FIXATION

TUBE TIDY CLIPS




The Tube Tidy Clips are intended for keeping tubes and lines arranged and maintained during patient monitoring.

The product can clip to blankets, pillows, sheets and opens to allow tubing to run through. Simply lay tubes required over hook/loop and seal to hold tubes in place.

Organises cords and drapes, useful in all specialist units and wards.

Can retain drapes and linens without the need of pins or fixers, thus eliminating the chance of personal injury.

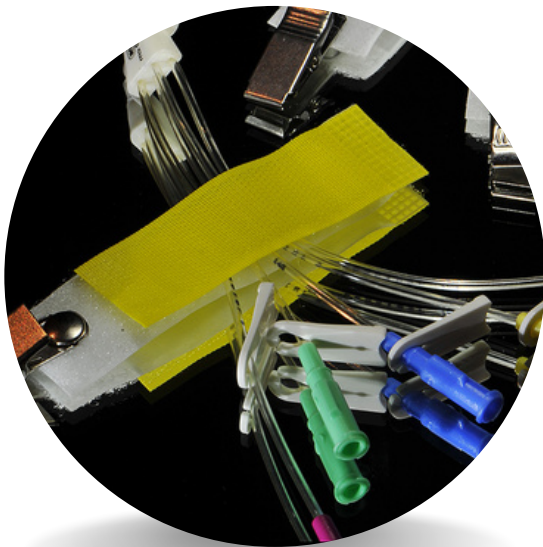

- Available in white or coloured for tube identification
- User is able to write on the product for identity markings
- Ease of use
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

The Tube Tidy Clips are used to organise and prevent damage to tubes and lines during patient monitoring. The coloured variants can be used for identification purposes to prevent incorrect tubes being applied by medical professionals.

Safety Information

- Do not apply overload weight to the clip. This increases the risk of the clip detaching from its position
- Attach the clip to only drapes and linen
- Discard the device should there be signs of deterioration
- Single patient use only

Medical Symbols



LATEX FREE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION

IMP CODE	DESCRIPTION	COLOUR	PACK QTY
TFP/800	Tube Tidy Clips	White	10
TFP/850.B	Tube Tidy Clips	Blue	20
TFP/850.R	Tube Tidy Clips	Red	20
TFP/850.Y	Tube Tidy Clips	Yellow	20

LARYNGECTOMY PROTECTORS

SOFSHIELD LARYNGECTOMY PROTECTOR




Sofshield Laryngectomy Protectors can be washed a maximum of 10 times, while still retaining filtration efficiency.


This offers a significant cost saving and longer product use over alternative protectors that use foam inserts, while also eliminating the need for costly additional cosmetic covers.

Sofshield is available in a variant of colours and sizes to suit all Laryngectomy and Tracheostomy patients.

Sofshield uses a silicone toggle system that eliminates the need to tie cords, therefore allowing patients with restrictive movement an easier fitting process.

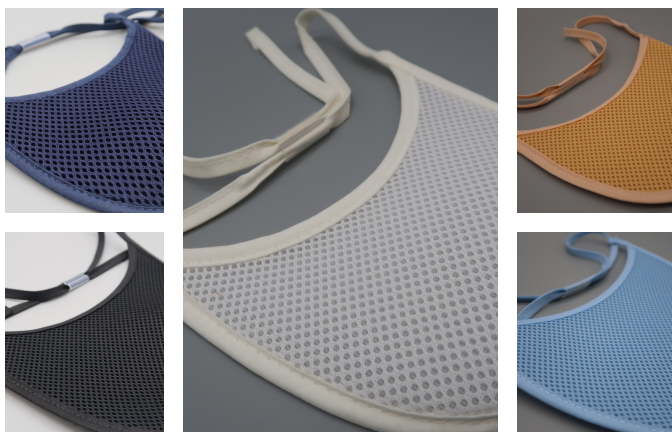

- Retaining filtration efficiency after washing
- Available on FP10 Prescription
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

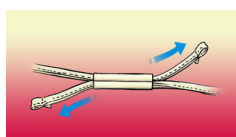

Intended Purpose

Sofshield Laryngectomy Protector is developed to help with the treatment for patients with a Laryngectomy. The protector is developed to filter, warm and humidify the incoming air.

Safety Information

- Do not use on broken or irritated skin
- Consult Instructions for use for washing information
- Single patient use only

Medical Symbols



LATEX FREE



Distributed by Insight Medical Products Ltd.

ORDERING INFORMATION



IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
AS/3669	Sofshield Protector -White	Large	291-9231	FTH006	10
AS/3764	Sofshield Protector -Blue	Large	305-5373	FTH008	10
AS/3765	Sofshield Protector -Beige	Large	305-5381	FTH007	10
AS/4051	Sofshield Protector -Black	Large	377-1722	FSM2621	10
AS/3698	Sofshield Protector -White	Small	291-9249	FTH009	10
AS/3777	Sofshield Protector -Blue	Small	305-5399	FTH100	10
AS/3778	Sofshield Protector -Beige	Small	305-5407	FTH010	10
AS/4063	Sofshield Protector -Black	Small	426-8199	FDf85026	10
AS/4056	Sofshield Protector -Light Blue	Medium	409-1526	DIRECT ONLY	10

LARYNGECTOMY PROTECTORS

CASCADE SHOWER PROTECTOR



Cascade is a soft waterproof protector made from a polyurethane material with a knitted polyamide backing that will absorb moisture from surrounding wet skin.

The neckband is soft and elastic to form a comfortable yet firm seal to prevent water seeping below the protector. The polyamide backing helps to prevent adhesion to the skin when the bib is wet reducing any potential breathing problems.

For added comfort the Cascade Protector can be worn over a Sofshield Laryngectomy Protector, however, this is not essential due to the backing material.

The protector is aesthetically pleasing white in colour and supplied in pairs to allow drying in between uses.



- Lay flat contour
- Washable (hand wash only)
- Available on FP10 Prescription
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

Cascade Shower Protector is developed to help with the treatment for patients with a Laryngectomy. Cascade is a waterproof protector to be worn around the neck. The material is soft and contoured to repel water spray away from a wide area around a stoma thereby improving patient confidence when showering.

Safety Information

- Replace at least every 6 months, or immediately if any signs of deterioration or discolouration in the material are observed
- Cascade may not be suitable for all types of deep stoma, if in doubt please consult your physician or nurse specialist
- Consult Instructions for use for washing information
- Single patient use only

Medical Symbols



LATEX FREE



Distributed by Insight Medical Products Ltd.

ORDERING INFORMATION



IMP CODE	DESCRIPTION	COLOUR	DRUG TARIFF	NHS CODE	PACK QTY
AS/3835	Cascade Shower Protector	White	314-2601	FDF85025	10 (5X2)

LARYNGECTOMY CARE


PURIFOAM LARYNGECTOMY FILTERS



Purifoam uses a high quality medical grade foam that protects the stoma from dirt and dust entering the airway, while warming and humidifying, the incoming air.

The filter is adhered to the skin using a skin friendly hypoallergenic pressure sensitive adhesive strip that has rounded edges for additional comfort.

Manufactured from beige coloured foam to give a natural aesthetic look. Purifoam is developed to have moisture retaining qualities, therefore ensuring good levels of absorbency and humidification to incoming air.

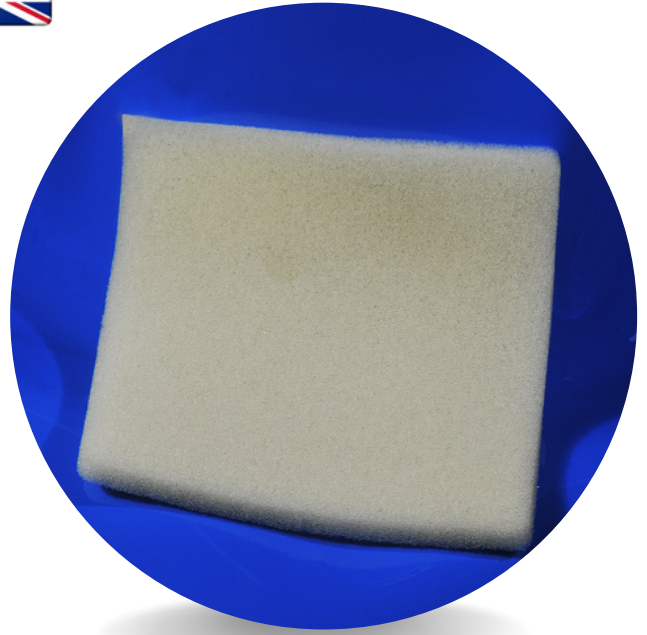

- Available on FP10 Prescription
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

Purifoam is developed to prevent dust and dirt entering the airway for patients with a Laryngectomy, while warming and humidifying, the incoming air.

Safety Information

- Change daily or more frequently if required
- Discard the device should there be signs of deterioration
- Single patient use only

Medical Symbols



LATEX FREE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION

NHS
Supply Chain

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/3853	Purifoam Filters	6.5x6.5cm	319-0386	FTH029	30

LARYNGECTOMY CARE

VOICE PROSTHESIS BRUSHES



Our Voice Prosthesis Brushes are manufactured with soft nylon bristles that prevent damaging the silicone inner surface of the prosthesis.

Brushes are tapered to provide effective cleaning of mucus and secretion build-up in the inner-lumen of the voice prosthesis.

The brush is fitted with a silicone safety sleeve to prevent over-insertion whilst in use.



- Silicone safety sleeve
- Polycarbonate handle for ease of use
- Available on FP10 Prescription
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

The Voice Prosthesis Brush is developed to ensure the cleaning of the tracheoesophageal voice prosthesis device fitted in certain laryngectomy patients as part of surgical voice restoration. The brush cleans the voice prosthesis and prevents damage which can lead to infections.

Safety Information

- The voice prosthesis may be cleaned while fitted, but extreme care must be used
- Do not push the brush right through the prosthesis
- Please seek help from your ENT Nurse or Speech Therapist should you be unsure of any element of using the device
- Bending and re-straightening the handle will cause undue stress, which can cause the handle to fracture
- We recommend using the brush a maximum of 5 times
- Do not use if allergic to silicone (Silicone Sheath Component)
- Discard the device should there be signs of deterioration
- Single patient use only

Medical Symbols



LATEX FREE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	PACK QTY
PB/50	Voice Prosthesis Brush	Standard	229-4395	10

BITE GUARD



Our Bite Guards are a protective shield used to protect teeth against accidental damage from contact with rigid instruments during surgical procedures.

Enhanced design to provide safer retention when fitted. Manufactured from medical grade TPE for pliability and patient comfort. The Bite Guard has smooth edges and compatible with most teeth configurations.

Coloured bright yellow to allow operator visibility. Can be supplied sterile or non sterile in an easy peel blister pouch.



- Shaped for optimum fit to the patient
- Radio-opaque, autoclavable on the 134°C cycle up to 3.15 minutes (non sterile version)
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

Bite Guard is used short term during investigative surgery to prevent the maxillary teeth becoming damaged from rigid instruments and devices handled by medical professionals during surgery.

Safety Information

- Do not use if the packaging is damaged
- Non-Sterile is autoclavable on the 134°C cycle up to 3.15 minutes
- Single patient use only

Medical Symbols



LATEX FREE

Non Sterile:



Sterile:



STERILE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION

DIRECT ONLY

IMP CODE	DESCRIPTION	SIZE	PACK QTY
BG/100	Non-Sterile Bite Guard	Adult	20
BG/200	Sterile Bite Guard	Adult	20

EYE SHIELDS



Eye Shields are used after surgery to protect the eye and to hold dressings in place while holes in the shield provide breathability and to prevent and reduce particles from entering and irritating the eye.

The Eye Shield is a single moulded component and has smooth edges for patient comfort. The raised dome on the Eye Shield prevents the device encountering the eye area which would be uncomfortable for the patient.

Supplied bulk for trade suppliers of Theatre Procedure Packs. We also supply individually sterile blister packs for hospitals and can use customers own artwork if required.



Intended Purpose

The Eye Shield is used for short term protection of the eye post ocular surgery. The shield is held in place with medical tape and will prevent the risk of injury from the patient inadvertently rubbing or accidentally bumping the eye post-surgery.

Safety Information

- Do not use if the packaging is damaged
- Single patient use only

Medical Symbols



LATEX FREE

Non Sterile:



Sterile:



STERILE



Manufactured by Insight Medical Products Ltd.

- Universal fitting available for dual fitting
- Available in Clear, Blue Tint or Opaque White
- Left and Right fitting in Opaque White only
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

ORDERING INFORMATION

IMP CODE	DESCRIPTION	PACKAGING	COLOUR	PACK QTY
ES/500/C	Universal Eye Shield	Bulk Non-Sterile	Clear	2000
ES/500/B	Universal Eye Shield	Bulk Non-Sterile	Blue Tint	2000
ES/500/W	Universal Eye Shield	Bulk Non-Sterile	White	2000
ES/500/WR	Right Fitting Eye Shield	Bulk Non-Sterile	White	2000
ES/500/WL	Left Fitting Eye Shield	Bulk Non-Sterile	White	2000
ES/100/C	Universal Eye Shield	Sterile	Clear	100
ES/100/B	Universal Eye Shield	Sterile	Blue Tint	100
ES/100/W	Universal Eye Shield	Sterile	White	100
ES/100/WR	Right Fitting Eye Shield	Sterile	White	100
ES/100/WL	Left Fitting Eye Shield	Sterile	White	100

AMNIOTIC PERFORATOR



The Amniotic Perforator is anatomically shaped to make the hook more versatile and easier to use. The perforator has been designed to make the membrane easily accessible by following the body's anatomical shape.

The single moulded design ensures optimum efficiency and additional prevention of injury to the patient when used.

The rounded tip and protected hook simply engages and ruptures the membrane, saving the need to prod or puncture, therefore reducing chance of injury to the foetus.



- Individually supplied in sterile peel pouch
- Manufactured in the United Kingdom
- Latex Free
- CE Marked

Intended Purpose

The Amniotic Perforator is used for the treatment of safely rupturing the amniotic membranes during the Amniotomy procedure. This is performed in transient use during vaginal examination and monitoring as part of labour management by a doctor or midwife.

Safety Information

- Do not use if the packaging is damaged
- Single patient use only

Medical Symbols



STERILE

LATEX FREE



Manufactured by Insight Medical Products Ltd.

ORDERING INFORMATION



IMP CODE	DESCRIPTION	SIZE	PACK QTY
KH/100	Amniotic Perforator	Standard	100

CAT BLADDER IRRIGATION CATHETERS

Manufactured as an OEM product for the veterinary trade.

Single use disposable catheter with a smooth closed top and two side eyes, all with an adjustable suture collar and welded luer connector.

Available in two materials, polypropylene (without stylet) and polyethylene (complete with introducing stylet) for longer term treatment. Available lengths of 130mm and 160mm.

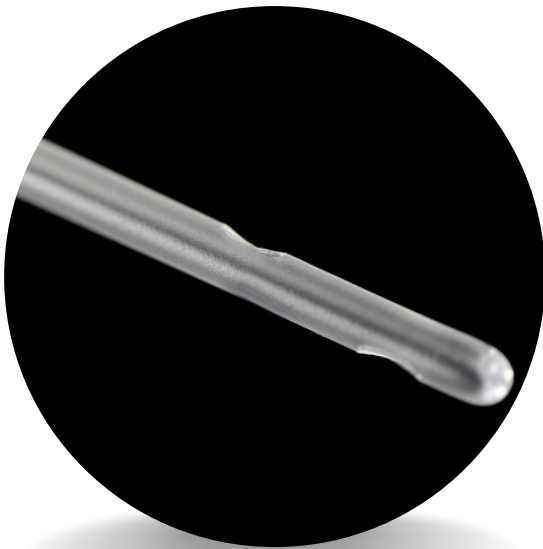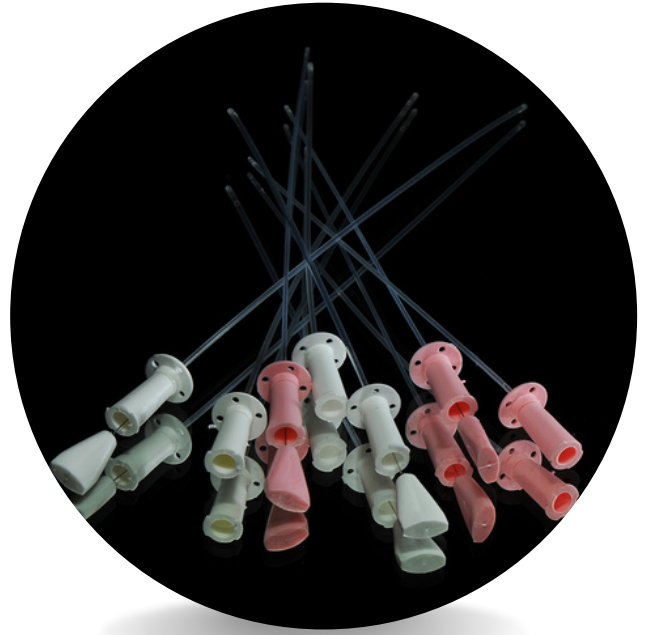

- Supplied packaged to customers requirements
- Blister packed, sterile with customers own artwork
- Bulk packed, unsterile for customers own packaging
- Manufactured in the United Kingdom
- Latex Free

ORDERING INFORMATION


IMP CODE	CATHETER SIZE (LENGTH)	MATERIAL	STYLET	LUER
CC/03/130P	3FG 1.0mm x 130mm	Polypropylene	No	Pink
CCS/03/130P	3FG 1.0mm x 130mm	Polyethylene	Yes	Pink
CC/03/160P	3FG 1.0mm x 160mm	Polypropylene	No	Pink
CCS/03/160P	3FG 1.0mm x 160mm	Polyethylene	Yes	Pink
CC/04/130W	4FG 1.3mm x 130mm	Polypropylene	No	White
CCS/04/130W	4FG 1.3mm x 130mm	Polyethylene	Yes	White
CC/04/160W	4FG 1.3mm x 160mm	Polypropylene	No	White
CCS/04/160W	4FG 1.3mm x 160mm	Polyethylene	Yes	White

DOG BLADDER IRRIGATION CATHETERS

Manufactured as an OEM product for the veterinary trade.

Single use disposable catheter with a smooth closed top and two side eyes and welded luer connector.

Available in polypropylene (without stylet). Available length of 500mm.

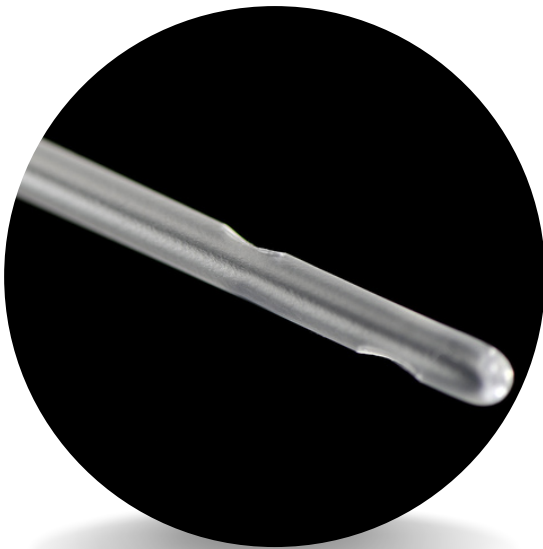

- Supplied packaged to customers requirements
- Blister packed, sterile with customers own artwork
- Bulk packed, unsterile for customers own packaging
- Manufactured in the United Kingdom
- Latex Free

ORDERING INFORMATION

IMP CODE	CATHETER SIZE (LENGTH)	MATERIAL	STYLET	LUER
DC/04/500/W	4FG 1.3mm x 500mm	Polypropylene	No	White
DC/06/500/Y	6FG 2.0mm x 500mm	Polypropylene	No	Yellow
DC/08/500/O	8FG 2.6mm x 500mm	Polypropylene	No	Orange
DC/10/500/R	10FG 3.3mm x 500mm	Polypropylene	No	Red


QUALITY & REGULATORY

Our non-sterile medical devices comply to the Medical Device Regulations (EU) 2017/745 and applicable harmonised standards. Our sterile devices meet the requirements to the Medical Device Directive 93/42/EEC, UK MDR 2002 and applicable harmonised standards under supervision of SGS. We are currently in transition for the sterile devices to be MDR Certified.

Our Quality Management System complies and meets the requirements of ISO 13485:2016 under supervision of SGS.

We are in progress of obtaining ISO 14001 -Environmental Management Systems certification under supervision of BAB.

For medical devices manufactured by Insight Medical Products, our EC Rep is as follows;

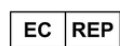

Advena Limited. Tower Business Centre, 2nd Flr,
Tower Street, Swatar, BKR 4013, Malta.

Reporting a Serious Incident under 2017/745 Medical Device Regulations

Contact Insight Medical Products Ltd. and the Competent Authority of the Member State in which the user and/or patient is established. This is for devices manufactured by Insight Medical Products which is clearly stated on each product page of this catalogue. If there is a complaint for devices distributed by Insight Medical Products, please refer to instructions with the device and contact Insight Medical Products for further actions.

PRESCRIPTION INFORMATION

Prescription Products are available from the following home delivery service providers:



Fittleworth Medical 0800 378 846



NWOS Homecare 0800 243 103



Patient Choice 0800 862 0686



ATOS Care 0800 783 1659


Also any pharmacies that offer a delivery service





INSIGHT
Medical Products Ltd


Unit 3, Priory Industrial Estate
Tetbury, Gloucestershire
GL8 8HZ
United Kingdom

 +44 (0)1666 500055

 info@insightmedical.net

Registered in England No: 349101

 InsightMedical

 @InsightMP

